


Brewer Task Lighting

LED Eco-Series

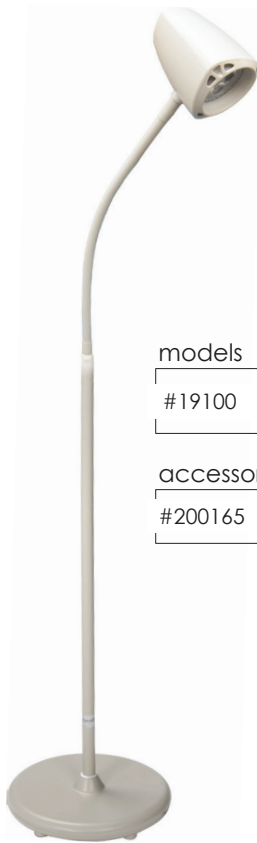

models

#19100 beige

accessories

#200165 five-caster base - beige

Light

- Standard one-year warranty
- Eco-friendly long-life LED bulb for lower cost of ownership
- Bright, best-in-class white light drives clinical efficiency
- Cool to the touch and easy to maneuver
- Head rotation: 170° Left/Right
- High intensity, low voltage, 12-Watt LED bulb included
- Brightness: 880 Lumens; 5" Beam dia. at 18" distance
- Color temperature: 4000° Kelvin
- Bulb life: Approx. 50,000 hrs

Arm

- Smooth, easy-to-clean flexible arm and head
- Flexible, "no-drift," 23" (58.4 cm) gooseneck
- Maximum Height: 66" (167.6 cm)

Base

- Power cord is detachable and exits the weighted base at the floor level
- Optional five-caster base available
- Base Diameter: 13" (33 cm)
- Cord length: 10' (304.8 cm)

Halogen 35 Task Series

Light

- Standard one-year warranty
- Crisp, white light
- Non-metallic head eliminates surface heat
- Dichroic film-coated reflectors filter out 80% of beam heat
- Heat sink technology removes bulb socket heat and keeps the housing cool
- High intensity, low voltage, 35-Watt bulb included
- Color temperature - 3000° Kelvin
- 1500 Footcandles @ 24" (61 cm)
- Bulb life: 4000 hrs

Arm

- Flexible, "no-drift," 24" (61 cm) gooseneck
- 540° rotation of the light head and gooseneck arm
- Maximum Height: 65.5" (166.4 cm)

Base

- Internal power cord exits the weighted base at floor level
- Optional five-caster base (Shown)
- Base Diameter: 13" (33 cm)
- Cord length: 9' 3" (281.9 cm)


models

#18100 beige

accessories

#001635 35-watt replacement bulb

#200165 five-caster base - beige

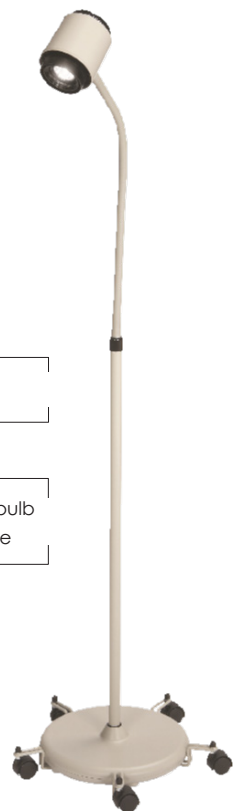

Brewer Task Lighting

Incandescent Task Series

Light

- Standard one-year warranty
- Crisp, white light
- Advanced reflector system provides shadow-free light
- Non-metallic head eliminates surface heat
- 25-Watt incandescent bulb included
- Bulb life: 500 hrs

Arm

- Flexible, "no-drift", 18" (45.7 cm) gooseneck offers superior beam alignment
- 540° rotation of the light head and gooseneck arm
- Maximum Height: 62" (157.5 cm)


Base

- Internal power cord exits the weighted base at floor level
- Optional five-caster base
- Base Diameter: 13" (33cm)
- Cord length: 9' 3" (281.9 cm)



models

#11500 chrome

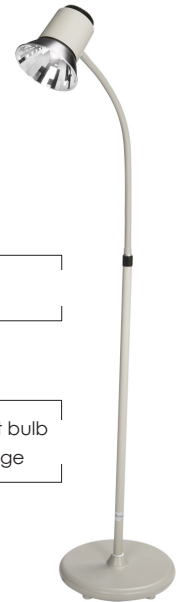

models

#17100 beige

accessories

#001525 25-watt replacement bulb
#200165 five-caster base - beige

Chrome Task Series

Light

- Standard one-year warranty
- 10 to 40-Watt incandescent bulb (not included)
- Some assembly is required

Arm

- Maintains its positioning properly over time
- Easily adjustable in any direction
- Durable chrome finish
- 11" (27.9 cm) gooseneck
- Maximum Height: 72" (182.9 cm)

Base

- Stable chrome plated base
- Base Diameter: 11" (27.9 cm)
- Grounded 3-wire plug
- Cord length: 12' (365.8 cm)

The Brewer Company is an ISO 13485:2003 certified company.

All Brewer products conform to US and Canadian electrical and medical industry codes and regulations and meet or exceed industry requirements.



All Brewer lighting comes with a 1-year warranty.

The Brewer Company reserves the right to improve or modify the design of its products consistent with new methods and techniques.

Call your local Brewer sales representative to learn about our complete product line.

The Brewer Company LLC
N88W13901 Main Street, Menomonee Falls, WI 53051
1.888.Brewer.1 www.brewercompany.com

Brewer